SUPERCRITICAL CPD MEMS DRYER FOR 6" (150mm) WAFERS AUTOMEGASAMDRI[®]-915B, SERIES B



- The Supercritical **Automegasamdri**[®]-915B, Series B System was developed by Dr. A. J. Tousimis and our team of long-standing associates within the tousimis[®] SAMDRI design group.
- Smooth operation and precise control are the Supercritical **Automegasamdri**[®]-915B's trademark; wafer holders and inserts are provided that allow anti-stiction processing up to 5 pieces of either 6", 4", 3", 2" wafers or 10mm square die per process run!

Tel: 301.881.2450 Fax: 301.881.5374 trc@tousimis.com • www.tousimis.com © 2022 @tousimis

®tousimis

SUPERCRITICAL Automegasamdri®-915B, Series B Catalog# 8785C (Process Various Wafer Sizes up to 6in - 150mm)

FEATURES

- Clean room static-free compatible design.
- Microprocessor controller allows for complete automatic processing.*
- All internal surfaces are inert to CO2 and ultrapure alcohols.*
- Repeatable operating parameters insuring reproducibility of results.*
- LED's instantly indicate process mode at a glance.*
- Integrated temperature controls.*
- Processes up to 5 x 6" (150mm) diam. wafers per run; comes with additional HF compatible wafer holders to process 5 x 4" (100mm), 5 x 3" (75mm),

5 x 2" (50mm) diameter wafers or 5 x 10mm square die (tousimis $^{\scriptscriptstyle B}$ HF compatible Wafer Holders*

can be used to transport and process your wafers and die from HF).

- All electronic components meet CE, UL and/or U.S. Military Specifications.
- Internal rupture disc built in for safety.*
- Static pressure control module pre-set for automatic safe pressure stability.
- An efficient adiabatic cooling system*, capable of cooling the chamber from ambient room temperature to near 0°C, in less than 3 minutes, maintains temperature automatically during both FILL and PURGE modes.*
- For added safety and convenience, U.S. Patented SOTER[™] Condenser* quietly captures exhaust and alcohols.
- Under-lit chamber with viewing window facilitates operator chamber status viewing.
- 0.08µm internal filtration system delivers clean filtered LCO2 to process chamber.*
- Unique chamber inserts* allow variance of chamber I.D. This maximizes efficiency in LCO₂ consumption and processing time.
- On-site Start-up and User(s) Training by a qualified tousimis® technician included.**


SPECIFICATIONS

- Cabinet: 14.25" (36.2cm) Width x 11.75" (29.8cm) Height x 25" (63.5cm) Length
- System Set-Up Area Required: 42" (107cm) Width x 30" (76cm) Depth.
- Chamber size: 6.50" I.D. x 1.25" Depth. Chamber Volume: 679 ml
- Temperature gauge range: -30°C to 60°C, Pressure gauge range: 0 to 2,000psi
- 120V/50-60Hz (Other voltage units also available. Please Inquire)
- LCO2 flow is controlled through Micro Metering Valves with Vernier handles for easily controlled flow rates.*

STANDARD ACCESSORIES

- LCO₂ High Pressure braided stainless steel inert Teflon[®] lined hose. 10ft (~3m) long for clean room operation (at a nominal charge, other lengths available upon request).
- Double T-Filter Assembly (#8785) pre-installed onto the chamber LCO₂ supply high-pressure hose. Flows LCO₂ twice thru 0.5 µm filters with 99.5% particulate retention prior to LCO₂ entering Automegasamdri[®]-915B.
- New SOTER[™] Condenser* connects to exhaust hose. Captures all alcohol and deadens exhaust noise.*
- Static free Exhaust Tubing for all exhaust outlets.
- Internal stainless and nickel scintered filtration systems incorporated to protect lines, wafers, and valves down to 0.4μm.
 Spare Chamber O-ring (3), Chamber Lamp (1), and 5A Slo-Blow Fuses (2).
- 4 Chamber Inserts.* Enables original chamber ID reduction down to smaller chamber ID sizes for a range including: 6", 4", 3", 2", and 1.25" !
- 6" (150mm), 4" (100mm), 3" (75mm), 2" (50mm) diameter Wafer Holders and 10mm square die holders included. Holders are HF compatible and can hold up to 5 wafers or die each.
- Free lifetime technical support consultation by our scientific staff.
- 2 year warranty on all parts and labor.


Note: Actual delivered model or accessories may vary slightly, as advancements are being constantly applied...

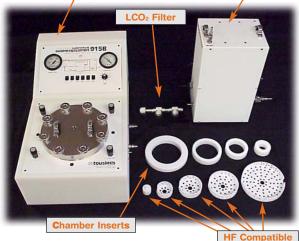

SOTER[™] Condenser

Wafer Holders

Automegasamdri®-915B

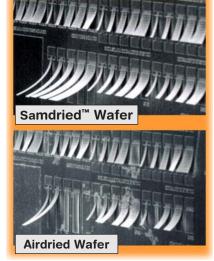

^{*} Most Automegasamdri[®] feature U.S. patents (# 6,493,964, #6,678,968) or patents pending.
** Available for most countries. Contact tousimis[®] for details.