


**OPTIONAL SOTER™ CONDENSER (#8777)**



Note: Overall System Set-Up Area Requirement with optional SOTER™ Condenser is 36”(W) x 30”(D) 92cm(W) x 76cm(D)

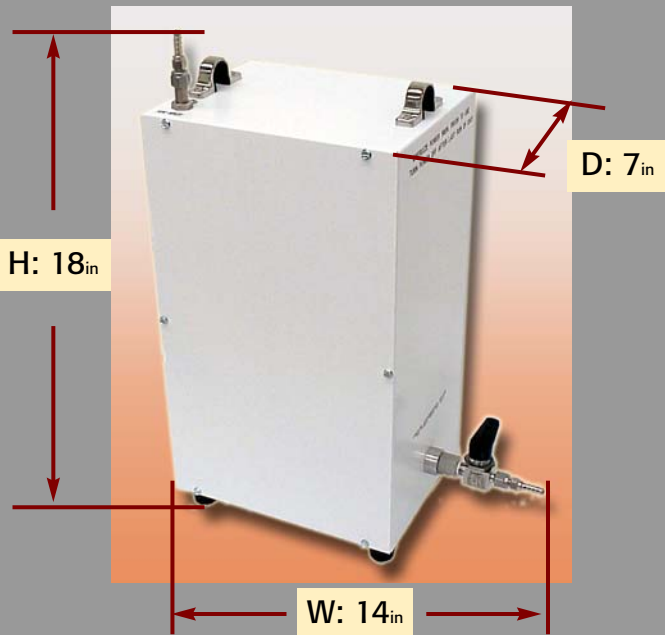

**8777 SOTER™ Condenser**

SOTER™ Condenser capable of quietly separating alcohol(s) from CO<sub>2</sub> exhaust discharge of your SAMDRI® chamber. Standard 120V grounded power source required. Comes complete with 10' static free 1/4" ID exhaust line, insulated exhaust connect from SAMDRI®, and clear drain line tubing.



**Optional SOTER™ Condenser Shown Connected to SAMDRI®-795**